

College of American Pathologists (CAP) GH5 Survey Data:

(updated 1/24)

The American Diabetes Association (ADA) recommends that “The A1C test should be performed using a method that is certified by the NGSP”. All laboratories performing HbA1c testing should participate in a fresh sample proficiency testing survey such as the College of American Pathologists (CAP). CAP GH5 data for the **third** survey of 2023 are summarized below. The NGSP target or reference values are based on replicate analyses using eight NGSP certified secondary reference methods.

Commentary by R. Little, Ph.D., NGSP Network Coordinator for the NGSP Steering Committee

Beginning in 2015 there are two CAP programs for HbA1c proficiency testing using fresh whole blood samples - GH2 and GH5. GH2 samples are shipped twice a year with three samples in each mailing as before. GH5 are shipped three times a year with five samples in each mailing. The three samples in each of the two GH2 mailings are also included in two of the GH5 mailings. Therefore, the NGSP follows the three GH5 surveys, which include all the samples used for both surveys.

In 2023, based on data from the GH5-C survey:

- **Bias from the NGSP target and variability ($\pm 2SD$) are shown in Table 1. Figure 1 shows the mean of each method compared to the NGSP/DCCT target (dotted line) for each sample. Symbols represent the mean of laboratory results for each method; error bars are $\pm 2SD$. The shaded rectangles reflect the current CAP acceptance limit of ± 6 . Graphs are in order by HbA1c level.**
- **Bias:**
 - The absolute mean bias for each method group ranged from 0 to 0.34% HbA1c. Two methods, Abbott Architect c and Vitros Chemistry, each had a mean bias $>0.3\%$ HbA1c for sample GH5-15.
- **Variability (CV)**
 - Method-specific, between-laboratory CVs ranged from 0.8% to 4.9%.
 - Unfortunately, there were 5 methods with CVs over 3.0% for one or more samples: Abbott Afinion AS100, Beckman AU HbA1c, Beckman AU Systems-Beckman reagent, Beckman UniCel Dx C Synchron, and Siemens DCA Vantage.
 - On the flip side, there were 9 methods with CVs $\leq 2\%$ for 5/5 or 3/3 samples– the Abbott Alinity ci and Architect c, Arkray Adams HA-8180 and HA-8190, Bio-Rad D-100, Roche Cobas c513, Sebia Capillarys 3, and Tosoh G8 and G11. 30% of participating laboratories are using these methods with CVs $\leq 2\%$.
 - Approximately 84% of laboratories are using methods with between-lab CVs $\leq 3\%$ for 5/5 HbA1c levels.
- **Pass Rates**
 - The current pass limit for the GH5 survey is $\pm 6\%$. The overall pass rates in this survey were 97.8%, 97.6%, 96.8%, 97.2% and 96.5% for samples GH-11 through GH-15. For individual methods, the lowest pass rate was 84.1% and the highest was 100%.
 - As expected, methods with small bias and low CVs will have the highest pass rates and, conversely, methods with large bias and/or high CVs will have the lowest pass rates.
- **The overall (all-method, all-lab) CVs for the last 33 surveys are shown in Table 2. Overall CVs for the 5 samples in the current survey ranged from 2.8% to 2.9%.**
- **Laboratories should avoid using methods with high CVs and/or consistent large bias.**

- **New Guidelines and Recommendations for Laboratory Analysis in the Diagnosis and Management of Diabetes Mellitus are now available: Diabetes Care, <https://doi.org/10.2337/dci23-0048>. Relevant to this CAP survey commentary: “The goals for imprecision for HbA1c measurement are intralaboratory CV <1.5% and interlaboratory CV <2.5% and ideally no measurable bias.”**

NOTE: The NGSP certification evaluates agreement of each method at the manufacturing site using one lot of reagents and calibrators, one instrument, and one application under optimal conditions. CAP precision reflects between-laboratory reproducibility, often with more than one lot of reagents and calibrators, and sometimes with different instruments (e.g. Siemens Advia instruments) and/or different applications (e.g. with or without sample pretreatment). In addition, if changes were made in the method just prior to NGSP certification, it is possible that not all participating laboratories in the field would have made the change at the time of the CAP survey. For these reasons, it is important that laboratories review not only the certification status of HbA1c methods but also their performance in the CAP survey over time (a good indication of field performance) when selecting or evaluating HbA1c assay methods.

TABLE 1: 2023 GH5-C

		GH5-11			GH5-12			GH5-13			GH5-14			GH5-15		
NGSP %HbA1c Reference Value (95% CI)		5.53 (5.49 - 5.57)			5.26 (5.22 - 5.31)			6.49 (6.44 - 6.53)			8.37 (8.33 - 8.42)			10.11 (10.06 - 10.16)		
Method name	no. labs (GH1-GH3)	Mean %HbA1c	Mean bias	% CV	Mean %HbA1c	Mean bias	% CV	Mean %HbA1c	Mean bias	% CV	Mean %HbA1c	Mean bias	% CV	Mean %HbA1c	Mean bias	% CV
Abbott Afinion 2	139	5.44	-0.09	1.9	5.26	0.00	2.0	6.44	-0.05	2.5	8.27	-0.10	2.2	10.08	-0.03	2.0
Abbott Afinion AS100	67	5.40	-0.13	2.6	5.25	-0.01	3.1	6.40	-0.09	2.4	8.31	-0.06	2.5	10.10	-0.01	1.1
Abbott Alinity ci series	189	5.46	-0.07	1.6	5.19	-0.07	1.9	6.39	-0.1	1.6	8.19	-0.18	1.5	9.81	-0.30	1.4
Abbott Architect c System	177	5.45	-0.08	1.7	5.20	-0.06	1.7	6.37	-0.12	1.5	8.18	-0.19	1.2	9.80	-0.31	1.2
ARKRAY Adams HA-8180 series	39	5.56	0.03	1.2	5.29	0.03	1.3	6.53	0.04	0.9	8.36	-0.01	0.8	10.14	0.03	1.0
ARKRAY Adams HA-8190 series	29	5.59	0.06	1.6	5.32	0.06	1.7	6.55	0.06	1.2	-	-	-	-	-	-
Beckman AU HbA1c Advanced	64	5.58	0.05	3.1	5.29	0.03	3.9	6.59	0.1	4.1	8.37	0.00	3.1	10.14	0.03	2.9
Beckman AU Systems - Beckman reagent	54	5.60	0.07	3.3	5.31	0.05	3.4	6.56	0.07	3.4	8.35	-0.02	3.3	10.03	-0.08	3.1
Beckman UniCel DxC Synchron Systems	40	5.53	0.00	2.6	5.23	-0.03	3.1	6.44	-0.05	2.5	8.34	-0.03	1.8	10.13	0.02	2.7
Bio-Rad D-10	90	5.59	0.06	2.9	5.21	-0.05	2.4	6.60	0.11	2.5	8.56	0.19	2.0	10.37	0.26	2.4
Bio-Rad D-100	138	5.43	-0.10	1.8	5.06	-0.20	2.0	6.40	-0.09	1.6	8.32	-0.05	1.8	10.04	-0.07	1.7
Bio-Rad Variant II Turbo 2.0	95	5.53	0.00	2.4	5.14	-0.12	2.5	6.47	-0.02	2.0	8.44	0.07	1.8	10.16	0.05	1.9
Bio-Rad Variant II	11	5.60	0.07	2.1	5.25	-0.01	2.0	6.63	0.14	1.9	8.52	0.15	1.4	10.21	0.10	1.9
Roche cobas c311	18	5.40	-0.13	3.0	5.09	-0.17	2.6	6.40	-0.09	2.5	8.30	-0.07	2.3	10.09	-0.02	2.5
Roche cobas c500 series	452	5.54	0.01	2.8	5.27	0.01	2.8	6.53	0.04	2.4	8.40	0.03	2.2	10.10	-0.01	2.5
Roche cobas c513	68	5.60	0.07	1.7	5.32	0.06	1.9	6.56	0.07	2.0	8.37	0.00	2.0	10.11	0.00	2.0
Roche COBAS Integra 400	24	5.52	-0.01	1.8	5.28	0.02	2.6	6.47	-0.02	2.0	8.37	0.00	1.7	10.13	0.02	2.6
Sebia Capillars 2 Flex Piercing	43	5.49	-0.04	2.3	5.23	-0.03	2.7	6.48	-0.01	1.9	8.29	-0.08	1.4	9.90	-0.21	1.9
Sebia Capillars 3 (CAPI 3) Tera/Octa	45	5.50	-0.03	1.4	5.25	-0.01	1.9	6.47	-0.02	1.5	8.29	-0.08	1.2	9.88	-0.23	1.3
Siemens Atellica CH-Enzymatic reagent	152	5.42	-0.11	2.8	5.16	-0.10	2.6	6.40	-0.09	2.4	8.24	-0.13	1.9	9.92	-0.19	1.5
Siemens DCA Vantage	262	5.53	0.00	3.9	5.30	0.04	3.2	6.50	0.01	4.1	8.42	0.05	4.9	10.17	0.06	4.8
Siemens Dimension ExL	175	5.56	0.03	2.7	5.24	-0.02	2.5	6.55	0.06	2.4	8.30	-0.07	2.1	10.15	0.04	2.0
Siemens Dimension Vista	134	5.41	-0.12	2.1	5.12	-0.14	2.1	6.36	-0.13	1.9	7.94	-0.43	2.5	9.84	-0.27	1.7
Tosoh G8 Automated HPLC	252	5.68	0.15	1.5	5.37	0.11	1.7	6.67	0.18	1.4	8.61	0.24	1.3	10.33	0.22	1.5
Tosoh G11 Automated HPLC	11	5.56	0.03	1.2	5.26	0.00	1.0	6.52	0.03	1.0	-	-	-	-	-	-
Trinity Biotech Premier Hb9210 HPLC	68	5.55	0.02	2.1	5.35	0.09	2.2	6.50	0.01	2.0	8.45	0.08	1.8	10.33	0.22	2.0
Vitros Chemistry Systems (5,1 FS; 4600; 5600; XT 3400; XT 7600)	209	5.42	-0.11	2.4	5.13	-0.13	2.2	6.43	-0.06	2.6	8.26	-0.11	2.7	9.77	-0.34	2.9

Gray shading indicates bias > 0.3% HbA1c or CV > 3.0% Note: these are arbitrary limits chosen to highlight methods with the highest bias and CV.

Figure 1: Bias and Variability from the NGSP Target

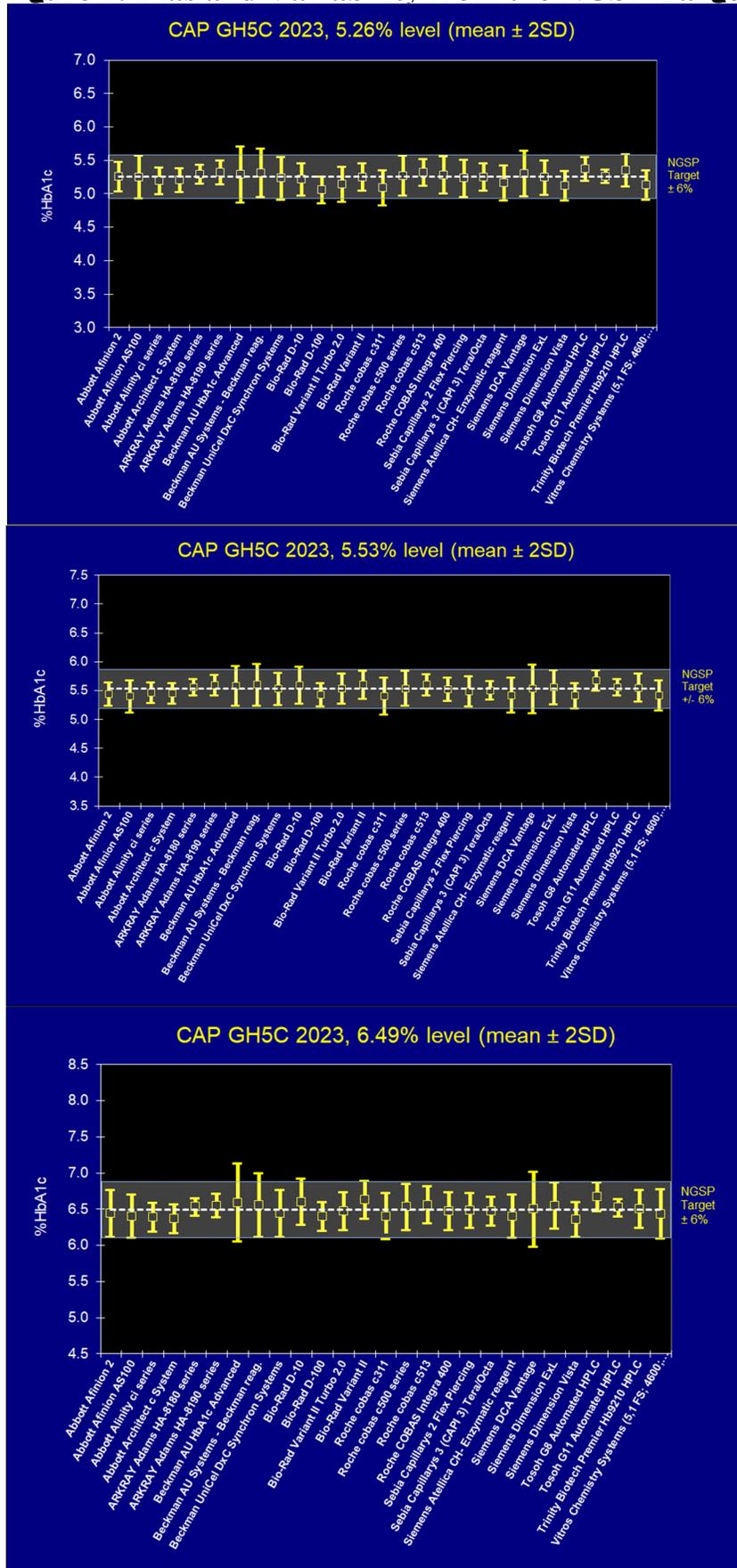


Table 2: Overall Variability for 2010-2023 for all GH participants

Mailing	Sample#	# of labs	Target	All method mean	S.D.	C.V.
A-2010	01	2573	5.9	6.03	0.23	3.9
	02	2566	9.8	9.73	0.39	4.0
	03	2581	7.4	7.43	0.31	4.2
B-2010	04	2693	5.2	5.34	0.21	4.0
	05	2691	8.7	8.67	0.33	3.8
	06	2685	6.3	6.37	0.23	3.5
A-2011	01	2652	8.5	8.58	0.28	3.2
	02	2645	5.4	5.52	0.20	3.5
	03	2649	6.4	6.51	0.21	3.2
B-2011	04	2877	6.3	6.36	0.24	3.8
	05	2872	7.6	7.69	0.29	3.8
	06	2871	9.2	9.28	0.34	3.7
A 2012	01	3298	5.6	5.62	0.20	3.5
	02	3316	9.4	9.44	0.37	3.9
	03	3301	7.2	7.28	0.29	3.9
B2012	04	3222	5.4	5.51	0.21	3.9
	05	3208	8.3	8.31	0.31	3.7
	06	3172	5.65	5.75	0.32	5.6
A 2013	01	2816	7.1	7.12	0.25	3.5
	02	2829	9.3	9.39	0.31	3.3
	03	2840	6.1	6.13	0.24	3.9
B2013	04	2912	8.1	8.04	0.31	3.8
	05	2907	5.3	5.33	0.20	3.8
	06	2908	6.1	6.17	0.24	3.9
A2014	01	3277	6.5	6.60	0.25	3.8
	02	3267	7.0	7.09	0.27	3.8
	03	3253	9.7	9.72	0.33	3.4
B2014	04	3278	6.58	6.64	0.23	3.5
	05	3273	8.39	8.45	0.30	3.6
	06	3266	5.65	5.67	0.21	3.6
A2015	01	3237	6.79	6.82	0.25	3.6
	02	3246	10.28	10.19	0.36	3.5
	03	3252	6.82	6.82	0.25	3.6
	04	2365	8.63	8.63	0.30	3.4
	05	2362	5.32	5.36	0.18	3.4
B2015	06	2379	5.84	5.87	0.2	3.5
	07	2392	11.71	11.68	0.44	3.8
	08	2402	9.53	9.5	0.33	3.5
	09	2386	5.04	5.08	0.17	3.4
	10	2403	7.38	7.35	0.26	3.5
C2015	11	3284	11.69	11.68	0.47	4.1
	12	3285	5.93	5.95	0.19	3.3
	13	3286	5.17	5.2	0.17	3.3
	14	2410	8.14	8.12	0.24	2.9
	15	2408	9.3	9.25	0.29	3.2
A2016	01	3358	5.32	5.33	0.16	3.1

	02	3365	9.17	9.21	0.28	3.0
	03	3357	5.31	5.33	0.16	3.1
	04	2425	12.03	12.12	0.40	3.3
	05	2419	5.94	5.96	0.16	2.8
B2016	06	2433	5.27	5.27	0.15	2.8
	07	2427	10.59	10.55	0.33	3.1
	08	2440	6.20	6.17	0.18	2.9
	09	2428	12.23	12.21	0.44	3.6
	10	2443	7.51	7.52	0.20	2.7
C2016	11	3377	9.11	9.08	0.24	2.6
	12	3402	6.01	5.99	0.16	2.6
	13	3372	11.71	11.69	0.39	3.4
	14	2432	5.02	5.02	0.16	3.2
	15	2442	7.58	7.58	0.19	2.5
A2017	01	3418	6.41	6.45	0.2	3.1
	02	3393	9.53	9.56	0.28	2.9
	03	3409	5.34	5.35	0.17	3.1
	04	2461	8.51	8.58	0.23	2.7
	05	2460	7.25	7.33	0.19	2.6
B2017	06	2446	7.42	7.49	0.22	3.0
	07	2445	5.2	5.22	0.19	3.6
	08	2450	8.31	8.33	0.21	2.5
	09	2443	10.33	10.38	0.31	3.0
	10	2457	5.87	5.91	0.19	3.2
C2017	11	3331	7.25	7.28	0.19	2.6
	12	3339	6.81	6.85	0.2	3.0
	13	3344	8.65	8.64	0.23	2.6
	14	2431	9.5	9.54	0.26	2.8
	15	2419	5.45	5.45	0.16	2.9
A2018	01	3371	7.15	7.16	0.21	3.0
	02	3343	5.19	5.20	0.19	3.6
	03	3369	8.42	8.39	0.24	2.9
	04	2466	9.79	9.75	0.29	2.9
	05	2473	6.12	6.13	0.18	3.0
B2018	06	2442	6.31	6.34	0.19	2.9
	07	2443	9.11	9.09	0.25	2.8
	08	2446	5.31	5.30	0.15	2.9
	09	2448	8.35	8.40	0.21	2.5
	10	2446	7.62	7.64	0.19	2.5
C2018	11	3311	9.37	9.31	0.27	3.0
	12	3291	4.92	4.96	0.16	3.2
	13	3300	6.04	6.09	0.18	3.0
	14	2463	5.28	5.32	0.17	3.1
	15	2469	8.08	8.11	0.24	3.0
A2019	01	3136	5.46	5.45	0.20	3.6
	02	3089	5.66	5.71	0.25	4.4
	03	3232	9.31	9.29	0.29	3.1
	04	2470	5.28	5.24	0.17	3.3
	05	2482	7.41	7.43	0.20	2.7

B2019	06	2462	6.41	6.44	0.18	2.8
	07	2460	8.60	8.66	0.29	3.3
	08	2461	5.42	5.41	0.16	2.9
	09	2467	7.38	7.44	0.19	2.6
	10	2457	9.75	9.77	0.31	3.2
C2019	11	3268	5.25	5.20	0.16	3.0
	12	3283	6.41	6.44	0.19	3.0
	13	3283	8.21	8.23	0.24	2.9
	14	2463	5.47	5.45	0.15	2.8
	15	2461	9.48	9.50	0.28	3.0
A2020	01	2862	7.97	8.06	0.22	2.7
	02	2976	5.89	5.92	0.16	2.7
	03	2962	5.13	5.11	0.15	2.9
	04	2311	7.40	7.43	0.22	2.9
	05	2328	9.17	9.17	0.25	2.7
B2020	06	2392	7.81	7.80	0.20	2.6
	07	2398	5.91	5.90	0.18	3.0
	08	2394	8.69	8.69	0.25	2.8
	09	2408	5.86	5.90	0.18	3.1
	10	2396	9.56	9.52	0.30	3.1
C2020	11	3197	5.48	5.46	0.16	2.9
	12	3204	8.26	8.26	0.22	2.7
	13	3209	5.15	5.11	0.14	2.8
	14	2403	10.05	10.07	0.28	2.8
	15	2411	5.99	5.97	0.16	2.7
B2021	06	2408	5.80	5.77	0.17	3.0
	07	2402	9.12	9.12	0.24	2.6
	08	2402	5.47	5.43	0.14	2.6
	09	2408	7.25	7.23	0.19	2.6
	10	2398	9.26	9.28	0.25	2.7
C2021	11	3156	5.21	5.16	0.14	2.7
	12	3176	6.68	6.66	0.16	2.4
	13	3159	---	8.32	0.34	4.0
	14	2382	5.98	5.97	0.14	2.4
	15	2392	10.12	10.15	0.26	2.6
A2022	01	3121	8.42	8.44	0.23	2.7
	02	3127	6.31	6.31	0.14	2.2
	03	3131	9.05	9.12	0.24	2.7
	04	2411	5.18	5.14	0.13	2.5
	05	2410	6.72	6.72	0.15	2.2
B2022	06	2393	5.93	5.91	0.15	2.5
	07	2393	8.92	8.96	0.22	2.5
	08	2405	5.44	5.39	0.14	2.5
	09	2390	7.33	7.34	0.19	2.6
	10	2382	9.37	9.42	0.23	2.4
C2022	11	2903	5.09	5.04	0.14	2.8
	12	2953	6.91	6.93	0.18	2.6
	13	2906	8.84	8.90	0.24	2.7
	14	2348	9.98	10.03	0.24	2.4
	15	2308	5.52	5.47	0.14	2.6
A2023	01	3012	7.58	7.51	0.20	2.7
	02	3007	8.88	8.84	0.22	2.5
	03	3012	5.82	5.78	0.17	2.9
	04	2374	10.00	9.97	0.23	2.4
	05	2382	5.24	5.17	0.15	3.0
B2023	06	2415	9.21	9.26	0.22	2.4
	07	2426	6.09	6.06	0.17	2.8
	08	2410	8.75	8.79	0.24	2.7
	09	2421	5.11	5.09	0.14	2.8
	10	2419	8.09	8.09	0.21	2.6

C2023	11	3089	5.53	5.52	0.16	2.9
	12	3084	5.26	5.24	0.15	2.9
	13	3088	6.49	6.49	0.18	2.8
	14	2413	8.37	8.33	0.23	2.8
	15	2420	10.11	10.04	0.30	2.9

CVs **below** 3.5% are highlighted in pink

CVs **below** 3.0% are highlighted in blue

HbAS samples is indicated in yellow

There is no assigned value for GH13 due to an unusual amount of variability found during value assignment. Interestingly, there was also a higher than usual variability seen in the all method CV%.